



MULTI FUNCTIONAL ORA DEVICE

User Manual

www.ora-redlighttherapy.com





Thank you for choosing our product.

To ensure proper usage, please read this user manual carefully before using the product.

Please keep it in a safe place for future reference.

If you transfer the product to someone else, please include this user manual.



User Manual

Table of Contents

SAFETY GUIDELINES	01
PRODUCT INTRODUCTION	02
PRODUCT STRUCTURE (Diagram)	03
PRODUCT PARAMETER (Data Table)	04
USAGE METHOD	05
USAGE METHOD OF COTTON PAD	05
MODE INTRODUCTION & USAGE METHOD	07
Mode 1: Cleaning	07
Mode 2: Eyes	07
Mode 3: RF	07
Mode 4: Cool	08
Mode 5: Daily	08
Mode 6: Weekly	08
DANGER WARNING	09
PROHIBITION ON WARNINGS	10
Prohibited Areas	10
Prohibited Users	10
INSTRUMENT CLEANING	11
INSTRUMENT STORAGE, MAINTENANCE, AND TESTING	12
TROUBLESHOOTING INSTRUCTIONS	13
GUARANTEE & DISCLAIMER STATEMENT	14



SAFETY GUIDELINES

1. Not suitable for individuals who have undergone facial plastic surgery or subdermal metal implants, as well as pregnant women.

2. Each mode has a specified duration of use. When a mode ends, the device will enter the preparation stage of the next mode. During the preparation stage, only the mode name will be displayed, and the white backlight will flash.

To officially enter the mode, press the mode button, and the blinking white backlight will become steady. After mode 6, the device will shut down. If there is no operation during the preparation stage, the device will automatically shut down after 1 minute, and the light will turn off.

Please follow the instructions and do not exceed the recommended usage time.

When using the device, hold the sensing decoration on the back shell and use the massaging head to make contact with the skin, forming a circuit. Please ensure that the skin is moist during use, otherwise the device will not work properly.

4. When using the device on someone else's face, it is necessary to form a circuit by also touching the user's body. Otherwise, the device may not function properly.

5. When using the CLEANING mode, combine it with makeup remover and toner. The massaging head is waterproof, but the body of the device is not waterproof.

6. When using this product, please keep your face moist or a cotton pad wet. Avoid dry rubbing.

7. This device is suitable for the face and neck. When using it, please avoid the larynx. Please do not use on the following areas: wounds, surgical sites, inside the mouth, eyelids, eyeballs, and the area around the eyes. This product should not be washed or soaked. After use, wipe the device head with a damp cotton pad.

8. Please store the device in a cool and dry place. Do not place it in the bathroom or other damp areas to avoid short-circuiting or other malfunctions.

9.

10. This product has a warming function. To avoid discomfort on the skin, it is important for individuals. Please use the product according to the instructions and do not exceed the recommended usage time.

12. Please do not attempt to modify or repair the device yourself.



Warning

The company shall not be held responsible for any personal injury or financial loss resulting from failure to follow the instructions in the user manual before using this product. Please ensure strict compliance with the instructions to avoid any potential risks.

PRODUCT INTRODUCTION

The Radio Frequency Beauty Instrument is an innovative home beauty device designed for daily facial and body treatments. It is suitable for daily care of both the face and body. This versatile product not only helps with skin whitening, skin softening, freckle reduction, and body shaping but also provides microcurrent lifting and skin tightening effects.

Additionally, it features a fast hot compress and cold compress function using radio frequency technology. This allows for better absorption of nutrients, resulting in hydrated, radiant, and elastic skin. Experience the outstanding benefits of this professional-grade beauty instrument for enhanced beauty and skincare.



PRODUCT STRUCTURE



 **Note**

In order to continuously improve product specifications, the company reserves the right to make enhancements. The actual product may vary from the description provided.

**PRODUCT PARAMETER**

Product Name	Multifunctional ORA Device
Product Model	W-618J
Net Weight	211 g
Rated Voltage	3.7 V
Voltage/Current	5 V / 2 A
Rated Power	10 W
Battery Capacity	2000 mAh 7.4 Wh
Strength Range	3 gear
Mode	6 modes

USAGE METHOD

1. Press and hold the power/gear switch button  for 1.5 seconds to start or shut down the device.
2. Press the mode switch button  to cycle through the 6 available modes.
3. Press the gear switch button  to cycle through 3 levels of intensity.
4. During operation, there will be no charging indicator. During non-operational charging, the 3-level gear light indicates the battery level. The gear light will flash in white to indicate the current charging level. Once the battery is fully charged, the white light will stay steady, and then the next gear light will flash in white to indicate ongoing charging.

This process continues until all three gear lights are steadily white, indicating that the device is fully charged. When the power is insufficient, the low-level light will blink at a frequency of 2Hz for 10 seconds as a reminder, after which the device will automatically power off and the white light will turn off.


5. Each time the device is powered on, it defaults to mode 1 with low intensity. After switching modes, the initial intensity remains low.

6. Press the button will produce a "beep" sound. There will be no indication during operation. After each mode's preset timing is completed, there will be two long "beep" sounds as a reminder.

USAGE METHOD OF COTTON PAD

1. Remove the Buckle-pressed ring of cotton pad from the top of the device.
2. Place a cotton pad on the treatment head and secure it using the Buckle-pressed ring of cotton pad. (If the Buckle-pressed ring of cotton pad cannot secure the cotton pad, please replace it with a thinner and larger cotton pad. Using a thick cotton pad may weaken the skin's sensation of warmth.)
3. Apply toner to moisten the cotton pad on the treatment head. (When using skincare products, ensure that they completely saturate the cotton pad on the treatment head. If the cotton pad becomes dry, please promptly replenish it with skincare products.)



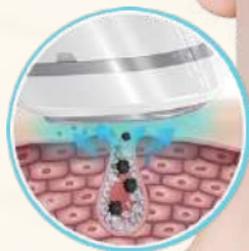

manual cleaning



Using Lumière's Device

***Deep clean and
purify pores!***

-  www.ora-redlighttherapy.com
-  Glendale, AZ USA
-  info@ora-redlighttherapy.com

MODE INTRODUCTION & USAGE METHOD

Mode 1: CLEANING

– Skin Cleaning and Nourishing Mode

- Working time: 4 minutes, with the head's LED light in orange color constantly on;
- The motor does not vibrate when there is no circuit, but vibrates constantly when there is a circuit;
- Adjustable motor vibration intensity with 3 levels;



Mode 2: EYES

– Eye Care Massage

- In this mode, the white backlight is constantly on, and it turns off when switching to the next mode;
- Working time: 3 minutes, head three-color light remains off;
- The motor operates in a 1-second on, 1-second off cycle (high-level intensity);
- The 3 gears indicate the intensity of EMS function; EMS function is available, operating in a cycle of 30ms on, is off, with an EMS frequency of 25kHz;
- No heating and cooling function; No RF heating function.



Mode 3: RF

– Skin Tightening and Wrinkle Reduction

- In this mode, the white backlight is constantly on, and it turns off when switching to the next mode;
- Working time: 4 minutes, with the head's LED light in red color constantly on;
- The motor is activated (high-level intensity);
- RF heating function, with 3 levels of low, medium, and high intensity in a cycling pattern, operating at a frequency of 2MHz.



Mode 4 : COOL

– Blue Light Cold Compress

- In this mode, the white backlight is constantly on, and it turns off when switching to the next mode;
- Working time: 6 minutes, with the head's LED light in blue color constantly on;
- The motor does not operate, and there is no adjustable intensity function for this mode. The intensity button is not functional;
- No EMS function, but cooling function is available with a temperature range of 11–15°C;
- No RF heating function.



Mode 5 : DAILY

– Daily Care

- In this mode, the white backlight is constantly on, and it turns off when switching to the next mode;
- Working time: 6 minutes, with the head's LED light in blue color constantly on;
- The motor is activated (high-level intensity);
- There are 3 levels indicating the EMS (Electrical Muscle Stimulation) strength;



- EMS function is available, alternating between strong and weak modes;
- No heating and cooling function, No RF heating function;

Mode 6 : WEEKLY

– Weekly Care

- In this mode, the white backlight remains on, and it turns off when switching to the next mode.
- Working time: 10 minutes, with the head's LED light in orange color constantly on;
- The motor does not operate.
- There are 3 levels indicating the EMS (Electrical Muscle Stimulation) and RF (Radio Frequency) heating intensity.
- EMS function is available, while there is no cooling function. RF heating function is available.



After the timed duration of each of the 6 modes, the current mode will end, and enter the preparation stage for the next mode (Mode 6 will directly shut down after the time duration). During this stage, only the mode name will be displayed, and the white backlight will flash. The mode will be officially entered by pressing the mode button, and the flashing white backlight will then become constantly on.

If no operation is performed during the preparation stage, the device will automatically power off after 1 minute, and the lights will turn off.



Danger Warning

1. Do not use near flammable substances such as alcohol, diluents, gasoline, nail polish remover, or styling sprays, as it may cause fire, electric shock, or bodily injury.
2. Do not use the power plug in the bathroom or in wet areas.
3. Do not let the device or charging interface get wet or clean with water, as it can result in electric shock, short circuit, or fire.
4. If there is an unusual odor or excessive heat, please stop using the device, as it may cause damage and lead to a fire.

5. Do not attempt to disassemble, repair, or modify the device yourself, as it may cause fire, electric shock, or bodily injury. If repairs are needed, please contact the after-sales service center.

6. When not in use or during the cleaning process, please make sure to turn off the power to avoid the risk of electric shock, short circuit, or fire.

7. Please disconnect the charging cable when the device is not in use. Do not remove or insert the USB charging head with wet hands, as it may cause electric shock, short circuit, or fire.

8. Keep the device away from heat source, direct sunlight, bathrooms, and other high-temperature or humid locations.

9. Avoid using the device in the same position for an extended period of time. It should be moved along the skin during use.

10. This product contains rechargeable batteries. Do not place the product in fire or heat it. It is also prohibited to charge, use, or store it in high-temperature environments, as it may result in overheating, fire, or damage.

11. When removing the USB charger, grasp the front end and avoid removing or inserting it with wet hands, as it may result in electric shock, short circuit, or fire.

12. Avoid placing it near heat-generating appliances, as it may cause overheating, fire, or damage. In case of any abnormal situation, immediately unplug the power plug. Do not continue using it if there are operational abnormalities, smoke, strange noises, or unusual smells.

13. Do not damage, modify, forcefully bend, pull, or twist the charging cable. Improper handling may cause the charging cable to break, leading to electric shock or fire.

14. This product must be used with the dedicated charging cable. Do not extend the charging cable on your own, as it may result in electric shock, short circuit, or fire.

15. Do not use it for purposes other than those described in the user manual. Otherwise, it may cause bodily injuries or product malfunctions

PROHIBITION ON WARNINGS

Prohibited Areas:

1. Do not use on wounds, surgical areas, oral cavity, genital area, mucous membranes, and other sensitive areas.
2. Do not use on skin with acne caused by inflammation.
3. Do not use on areas of the skin that are inflamed and suppurating.
4. Do not use on areas of the skin that are itchy, swollen, or painful due to skin cancer.
5. Do not use on eyelids, eyeballs, or areas around the eyes (except for eye care mode).

Prohibited Users:

1. Individuals with implanted medical devices such as artificial heart-lung machines, cardiac pacemakers, and electrocardiographs.
2. Individuals who are pregnant, menstruating, or breastfeeding.
3. Individuals taking sedative medications.
4. Individuals undergoing treatment.
5. Individuals with pathological spots caused by physical stimulation.

6. Individuals with dilated capillaries due to long-term use of corticosteroid hormones or liver dysfunction.
7. Individuals with atopic dermatitis.
8. Individuals with tuberculosis, fever, hypertension, or infectious diseases.
9. Individuals with congestion, blood disorders, visceral diseases, or allergies.
10. Individuals with acute illnesses during the growth process.
11. Individuals with impaired mental faculties or disabilities.

- Skin Tightening
- Freckle Removal
- Oil Control
- Anti-aging

INSTRUMENT CLEANING

After each use, please clean and maintain the instrument promptly to ensure cleanliness.

1. Before cleaning, make sure the power is turned off. Avoid using wet towels or paper towels for maintenance while the power is on, as this may cause irritation.
2. Do not immerse the product in water or wash it with water, as it may damage the instrument. Instead, use a dry cloth or a soft cloth dampened with a small amount of soapy liquid to wipe it.
3. Clean the Buckle-pressed ring of cotton pad. If there is dirt on the ring, rinse it directly with water.
4. Do not use corrosive solutions (such as alcohol, hydrogen peroxide, gasoline, or nail polish remover) to wipe the device. These substances may cause product failure, damage, or discoloration.



INSTRUMENT STORAGE, MAINTENANCE, AND TESTING

1. When not in use for a long period of time, fully charge the instrument and clean it with a clean cloth or cotton pad before storing.
2. Do not place the product near heating devices, in bathrooms, or in high-temperature and humid areas. Avoid direct sunlight or unstable places. Store the
3. product in a place where children cannot reach it.
If the product has not been used for a long time and
4. needs to be used again, make sure to check and confirm that the product is in a normal and safe operating condition before use.
If the product exhibits any abnormalities, please refer to the troubleshooting instructions of the user manual for preliminary checks. If the product still cannot be used
5. after preliminary checks, please contact us promptly.



Deep Clean



EMS



Tight Skin



Eye Care



Red Light
Therapy + RF



Cold Therapy
Blue Light



TROUBLESHOOTING INSTRUCTIONS

SYMPTOM	POSSIBLE CAUSE	SOLUTION
Unable to power on	No long-press of the power button	Press and hold the power button for at least 1.5 sec
	Insufficient battery charge	Please charge it in time, and use it after it is fully charged
Ionic strength too high	Setting level too high	Adjust the setting level using the mode button
Low ionic strength (No sensation)	Setting level too low	Adjust the setting level using the mode button
	Insufficient skin moisture	Increase skin hydration by using skincare products
	Makeup cotton too thick	Replace with a thinner makeup cotton
Overheating of thermal sensation heating	Position for a long time	Same area for more than 3 seconds
	Slow sliding speed	Use it slightly faster during operation
Heating of the device body	Using continuously multiple times	Turnoff the device and stop using it,then let it sit for at least 1 hour
Unable to vibrate	COOL - Blue Light Cold Compress DAILY - Daily Care	These two modes do not have vibration

GUARANTEE & DISCLAIMER STATEMENT

Thank you for choosing ORA Red Light Therapy. Before proceeding with the use of our services, please carefully read and understand the following Guarantee & Disclaimer statement.

1. Guarantee & Warranty: We at ORA Red Light Therapy are committed to providing high-quality products and services to our valued customers. We guarantee that our products have been manufactured with the utmost care and attention to detail. In the unlikely event that you encounter any issues with our products within 30 days of purchase, please contact us immediately at info@ora-redlighttherapy.com, and we will work diligently to resolve the matter to your satisfaction.

Disclaimer of Warranties: (a) The use of our services, products, and any related website is at your sole risk. Our services, products, and website are provided on an "as is" and "as available" basis. We make no representations or warranties of any kind, express or implied, including but not limited to, warranties of merchantability, merchantable quality, fitness for a particular purpose, durability, title, and non-infringement. (b) We do not guarantee, represent, or warrant that your use of our services will be uninterrupted, timely, secure, or error-free. (c) We do not warrant that the results obtained from the use of our services will be accurate or reliable. (d) We do not warrant that our website or any third-party websites linked to our website are free from viruses or other harmful components. You acknowledge that the use of our website and any third-party websites is at your own risk.

Limitation of Liability: (a) ORA Red Light Therapy, its directors, officers, employees, affiliates, agents, contractors, interns, suppliers, service providers, or licensors shall not be liable for any injury, loss, claim, or any direct, indirect, incidental, punitive, special, or consequential damages of any kind, including, without limitation, lost profits, lost revenue, lost savings, loss of data, replacement costs, or any similar damages, whether based in contract, tort (including negligence), strict liability, or

otherwise, arising from your use of any of our services or any products procured using our services, or for any other claim related in any way to your use of our services or any product, including, but not limited to, any errors or omissions in any content or any loss or damage of any kind incurred as a result of the use of our services or any content (or product) posted, transmitted, or otherwise made available via our services, even if advised of their possibility. (b) In no event shall our total liability to you for all claims related to our services exceed the amount paid by you, if any, for accessing or using our services during the twelve months prior to the claim or cause of action arising.

4. Compliance with Laws: You agree to comply with all applicable laws and regulations in connection with your use of our services and products.

5. Indemnification: You agree to indemnify, defend, and hold harmless ORA Red Light Therapy and its parent, subsidiaries, affiliates, partners, officers, directors, agents, contractors, licensors, service providers, subcontractors, suppliers, interns, and employees from any claim or demand, including reasonable attorneys' fees, made by any third-party due to or arising out of your breach of these Terms of Service or the documents they incorporate by reference, or your violation of any law or the rights of a third-party.

By accessing our website or using our services, you acknowledge that you have read, understood, and agreed to the terms and conditions set forth in this Guarantee & Disclaimer statement. If you do not agree with any part of this statement, please refrain from using our services and exit our website immediately.

For any questions or concerns regarding this Guarantee & Disclaimer statement, please contact us at info@ora-redlighttherapy.com.

Thank you for your understanding and cooperation.



Are you ready to take your wellness journey to the next level?

Connect with us today to explore our range of products and discover the transformative potential of red light therapy.



Contact Us:

 www.ora-redlighttherapy.com

 Glendale, AZ USA

 info@ora-redlighttherapy.com

 Ora Red Light Therapy

Copyright © 2025