

## Operation

<b>ON / OFF</b>	Long press 3 Seconds to turn on/off
<b>Automatic mode</b>	Upon startup, the device enters automatic mode and activates a comprehensive 20-minute light therapy cycle. The cycle consists of three stages: Stage 1 (10 mins) - Infrared (invisible) and red-light therapy for improved circulation and collagen production Stage 2 (5 mins) - Blue light therapy for acne reduction and skin purification Stage 3 (5 mins) - Orange light therapy for increased skin radiance and overall skin health This automated feature allows users to effortlessly enjoy a complete light therapy session in the comfort of their own home.
<b>Manual mode</b>	While in operation, simply select the color and press the "T" key to adjust the working time of the device. With each short press of the key, the working time changes between three options: 15 minutes, 20 minutes, and 30 minutes. This convenient feature allows users to easily customize their light therapy sessions based on their individual needs and preferences.
<b>Adjust Intensity</b>	In working condition, short press "+" "-" key to adjust the brightness (intensity) of the lights.
<b>Charging</b>	When the battery indicator on your device starts blinking, it means that the battery is running low and needs to be charged. Please connect your device to a power source as soon as possible to ensure uninterrupted use.
<b>Usage</b>	It is recommended to use skin rejuvenation mask 3-5 times a week

## People unsuitable to use

1. This product has no fireproof function, please keep away from fire or corrosive liquid;
2. Do not disassemble or repair this product by yourself;
3. Please keep it at a place out of reach of children;
4. This product has waterproof function;
5. Do not press it heavily so as to avoid deformation and failure.

## Troubleshooting

Indication	Solution
<b>Device does not turn on</b>	*Long press 3 Seconds on ON/OFF switch to turn on *The device needs to be charged
<b>Battery low voltage</b>	The device needs to be charged
<b>USB plug/socket is bad</b>	*Re-insert *Replace the USB plug/socket
<b>Auto power-off</b>	*Long press 3 Seconds on ON/OFF switch to turn on *The device needs to be charged

# User Manual



# Beauty LED Mask

## Please Read This Entire Manual Before Using

Thank you for purchasing the spectrum skin rejuvenation mask, please read this manual carefully before using the mask. Please pay special attention to the safety matters so as to operate and use the mask correctly.

This manual provides users with safe and effective ways of operation. After reading, please keep it properly for future reference.

## Technology Introduction

This product uses the light guide material technology of nanotechnology to form different nano unit light waves.

Blue light: 415nm


Yellow light: 605nm

Red light: 630nm

NIR: 850nm

## Product List


- ① The mask x1
- ② Straps x2
- ③ Controller x1
- ④ Cable x1
- ⑤ USB charging cable x1
- ⑥ Eye protector x2
- ⑦ User manual x1

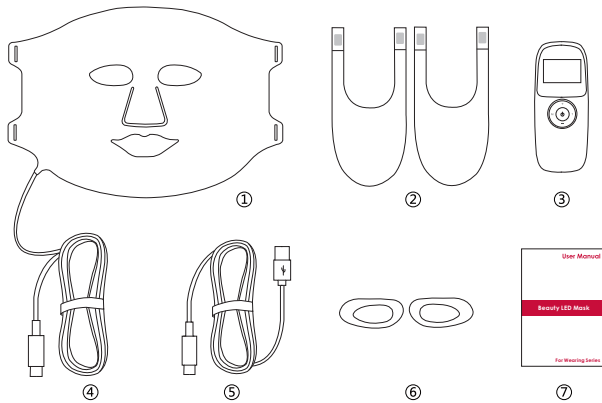

## Product Specification

<b>Product Name</b>	Beauty LED Mask
<b>Working voltage</b>	3.7V
<b>Charging voltage</b>	DC5V
<b>Rated power</b>	15W
<b>LED</b>	92pcs 4-chip
<b>Product Size</b>	315*219*5mm/12.4" x 8.6" x 0.2"
<b>Net Weight</b>	≈230g/0.5lbs
<b>Wavelength</b>	Blue light: 415±5nm,
	Orange light: 605±5nm
	Red light: 630±5nm,
	NIR: 850±5nm,
<b>Cell Rated Capacity</b>	2000mAh
<b>Material</b>	Food Grade Silicone


## Controller

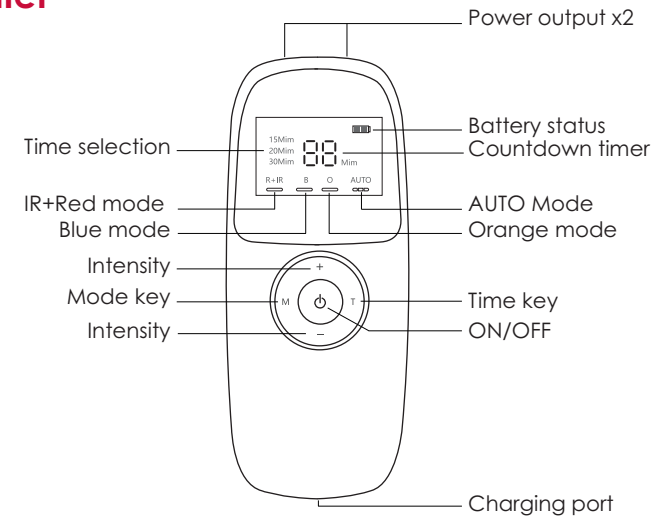

## Quick Start Guide



**Step 1:**  
Remove makeup and clean your face.  
Apply serum if desired.

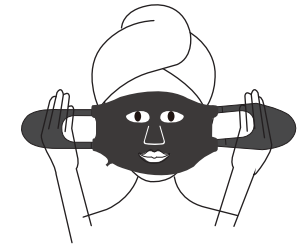

**Step 2:**  
Secure the device to your head using the velcro straps.



**Step 3:** Long press the power button for 3s to start the device. It runs automatically for 20 minutes on AUTO mode.



**Step 4:**  
Remove the mask and apply moisturizer to the skin.